Bedienungsanleitung
OPERATING MANUAL
NOTICE D'EMPLOI
BEDIENINGSHANDLEIDING
MANUAL DE OPERACIÓN

Elektronische Säuglings- und Kleinkindwaage
Electronic baby & toddler scale
Pèse-bébé/bambin électronique
Báscula electrónica de bebés y niños pequeños
Elektronische weegschaal voor baby's & peuters

M112800

M112800-160316-Rev002-UM-de-e-es-fr-nl

ADE GmbH & Co.
Hammer Steindamm 27 – 29

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☎ +49 (0) 40 432776-0
☎ +49 (0) 40 432776-10
✉ info@ade-hamburg.de

DE	Einleitung	3
	Allgemeine Sicherheitshinweise	3
	Beschreibung	4
	Waage aufstellen	6
	Spannungsversorgung	6
	Bedienungshinweise	6
	Reinigung und Pflege	7
	CE-Konformität	9
	Technische Daten	9
EN	Introduction	11
	Basic safety precautions	11
	Description	12
	Set up of the scale	14
	Power supply	14
	Operation	14
	Cleaning & maintenance	15
	CE-conformity	17
	Technical data	17
FR	Introduction	19
	Précautions élémentaires de sécurité	19
	Description	20
	Préparation de la balance	22
	Alimentation	22
	Utilisation	22
	Nettoyage et entretien	23
	Conformité CE	25
	Caractéristiques techniques	25
ES	Introducción	27
	Precauciones de seguridad básicas	27
	Descripción	28
	Preparación de la báscula	30
	Alimentación eléctrica	30
	Funcionamiento	30
	Limpieza y mantenimiento	31
	Conformidad CE	33
	Datos técnicos	33
NL	Inleiding	35
	Elementaire veiligheidsmaatregelen	35
	Beschrijving	36
	Instellen van de weegschaal	38
	Voeding	38
	Bediening	38
	Reiniging & onderhoud	39
	CE-conformiteit	41
	Technische gegevens	41

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Fachhändler:

Dear ADE customer,

Thank you for choosing our ADE electronic baby & toddler weighing scale. You purchased a high quality product for weighing babies & toddler for the purpose of monitoring and diagnosis.

Content

Description: **Electronic baby and toddler scale**

Model: **M112800**

1 x weighing scale base

1 x weighing tray

1 x weighing plate for toddlers

2 x AA batteries

Operation manual

Basic safety precautions

- Please read this manual carefully to get familiar with the features and operation before use.
- The usage and maintenance of the scale should to be done only by trained personnel.
- Ensure that the baby lies on the tray as relaxed as possible and without many movements.
- Never leave the baby unattended while on the scale.
- This product is intended for infants and children (toddlers) only.
- Use the scale on a firm and flat surface.
- Make sure that only the supplied type of battery is used. – Please see the operating instructions.
- Move all weight off the scale before use.
- Do not disassemble the scale as incorrect handling may cause injury.
- Replace the batteries when „Lo“ is being displayed. This is an indication of low battery.

Description

Weighing scale base

Weighing tray

Weighing

toddler

plate for

Description of the control panel


Control panel:	Description:
ON / OFF	Turn on and zero setting – automatic switch off
HOLD	Activates the auto hold function – to unlock press again
TARE	Activates the tare function

Display symbol:	Description:
◀	Zero position
Hold	Auto hold
kg lb	Weighting units kg or lbs

Set up of the scale

The scale consists out of 3 parts – weighing base, baby weighing tray and toddler weighing plate.

The baby weighing tray or the toddler weighing plate is fixed and released via the turning knob underneath the baby weighing base. Please make sure, that the tray or plate is completely and securely fixed before use.

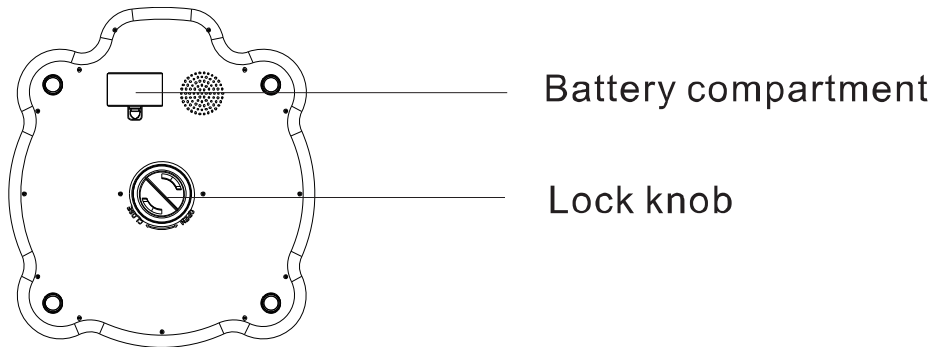Power supply

The scale is powered by 2 x AA alkaline batteries. Turn the scale around and open the battery compartment. Insert the enclosed AA alkaline batteries. Please observe the correct polarity as shown inside the battery compartment.

Replace the batteries when „Lo“ is being displayed. It indicates low battery power or use the mains connection.

Operation

Weighing

3. Start the scale by pressing the **ON/OFF** button. After a function test the last measured weight will be displayed for 3 seconds followed by 0.00 and ◀. The scale is ready for use and is set automatically to zero.
4. The baby can now be placed on the scale and the weight will be shown on the display in 10 g steps.

Auto hold function

Place the baby on the scale. The display will flush until the actual weight has been taken and **HOLD** will be shown on the display.

To cancel the hold function press the **HOLD** button on the key board.

The scale will switch of after approx. 60 seconds, when not in use.

Tare function

The tare function can be used to disregard an additional weight e.g. a nappy or towel on the weighing

tray surface. The tare function is operated as follow: Place the object on the tray surface and press the **TARE** button. The additional weight will be subtracted and 0.00 will be shown on the display. A baby can now be placed on the scale surface and the net weight will appear on the display. The maximum tare weight must not exceed 500g.

Unit switch function

Press the **ON/OFF** button twice during the function test. By pressing again the **ON/OFF** button you can select the requested units **kg** or **lbs**. The **kg** or **lb OZ** symbol will appear in the display to indicate the chosen weighing unit.

Recall of memorized weight

Start the scale via the **ON/OFF** button and the last measured weight will be displayed for approx. 3 sec. in order to determined weight changes of the infant.

Cleaning & maintenance

Cleaning instructions

Clean your scale with a soft moist towel or conventional disinfectant only. Do not use any aggressive cleaning solutions.

Warning signals and error messages

3. Why does the display read 'lo'?

Low battery. Replace the batteries.

4. The message displays 0-Ld

Overload warning; remove the weight immediately; otherwise, permanent damage to the scale will occur.

Safety instructions

The scale should never get wet. Since this can damage the electronics.

Waste disposal

Electronic devices are not domestic waste.

If the scale is out of use the responsible waste disposal association will be glad to advise you on the measures necessary to dispose of it properly.

Do not throw used batteries into domestic waste. Please inform yourselves at the responsible waste disposal association how to dispose used batteries.

Storage and transport

Keep all original parts and pieces for an eventual sending back of the scale in order to avoid any possible damage during transport which automatically leads to no warranty.

Disconnect all cables before transport to avoid unnecessary damages.

Warranty

A two year warranty period from date of purchase will be guaranteed for either repair or replacement once it is proven that the defects are caused by poor material or poor workmanship (please keep the receipt for presentation). All moveable parts such as batteries, cables, plug-in power supply units, rechargeable batteries etc. are excluded. Warranty does not apply to normal appearance of usage or sign of wear or any damages caused by accidents or misuse. Eventual warranty- or liability claims are subject to using the original ADE accessories or spare parts. Products which have been opened by unauthorized persons are automatically out of warranty. Customers abroad are requested to contact their local dealers in the event of any warranty becoming applicable.

CE-conformity

Declaration of conformity by the manufacturer

These scales have been manufactured in accordance with the harmonised European standards.
They conform to the regulations of the EC-directives listed below:

- EMC Directive 2004/108/EC regarding electromagnetic compatibility
 - Low Voltage Directive 2006/95/EC
 - RoHS Directive 2011/65/EC
- in the currently valid versions.

This declaration loses its validity if modifications are made to the scale without our approval.

Hamburg, August 2015

ADE GmbH & Co.

Hammer Steindamm 27 – 29
D-22089 Hamburg

Technical data

Capacity:	50 kg
Graduation:	10 g
Weighing unit selectable:	kg, lb
Battery operation	2pcs 1.5V AA alkaline batteries
Weight & Dimensions:	
Scale (LxWxH):	631 x 362 x 111 mm
Baby tray (LxWxH):	631 x 320 x 60 mm
Weight:	approx. 3,1 kg

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